Biofeedback & EMG-getriggerte Stimulation
Biofeedback & EMG-triggered stimulation

DuoBravo N

Medizinprodukt / Medical Device

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Progr.: PC4 - STIM-EMG-ETS konfigurierbar

Content

Please read the Operators Manual carefully
before using the device and pay attention
to all references !

In the following you'll find all information directly necessary for performing the treatment
with grey background.

	Page
I. Biofeedback- / ETS-device DuoBravo	
Delivery & Accessories / Medical Device Reporting	44
Technical Data	45
Explanation of Symbols	46
Safety Adjustments	46
Introduction – What is DuoBravo ?	47
Indication	48
Contraindication	49
Warnings	50
Advice on Application of Electrodes	51
Maintenance	52
Warranty	52
II. Display and Operating Elements	
Operating Elements	53
Battery Compartment	53
Connections	54
The Display	55
The Info-Line	56
III. Preparation of DuoBravo	
Basic Settings	57
Setup Real Time Clock	58
Setup Language	58
Preparation of Therapy	59

Content

	Page
IV. Therapy with DuoBravo	
The Different Phases	60
Use of pre-set programs	61
Therapy in a phase with pure stimulation	62
Therapy in an EMG-phase	63
Therapy in an ETS-phase	64
Therapy in combined programs	64
V. Additional Options for Configuration of DuoBravo	
Configuration of Free Programs	65
Parameters of Free Programs	66-69
Programming PC1	66
Programming PC2	67
Programming PC3	68
Programming PC4	69
Lock-Button	70
Threshold	70-71
Setting of threshold before start of treatment	70
Automatic adjustment of threshold during treatment	71
VI. Therapy-Control	
Therapy-Control	72-73
Statistics	74-75
Trouble Shooting Guide	76
VII. Programs	
Programs Neurology	78-82
Anatomy	83
Configuration	84
Waste Disposal / Recycling of Batteries / Literature / Copyright	85
Electromagnetic Compatibility	89-91
Graphics	92

Delivery & Accessories

Included in the delivery:

REF: 010E-090: 1 unit **DuoBravo N** (Neurology), 2 dual pin lead wires to connect the electrodes, 1 single pin reference lead wire, 8 self-adhesive electrodes (40 x 40 mm), 1 stand, 4 batteries type AA and 1 operators manual

Accessories:

REF

011E-006: 1 unit self adhesive electrodes 40x40 (4 pieces)

011E-007: 1 unit self adhesive electrodes 40x80 (4 pieces)

Spare parts:

REF

012E-120: 1 dual pin lead wire to connect the electrodes

012E-086: 1 carrying case

012E-083: 1 single pin reference lead wire

052F-805: 4 batteries type AA

012E-082: 1 operators manual **DuoBravo**

DuoBravo shall only be used with original accessories.

With regular average use, the expected service life of **DuoBravo** is 5 years. The battery life is approximately 12 hours or more at 20 mA per channel. The batteries can be stored for approximately 2 years. The service life (expiry date) of electrodes is printed on each package.

Medical Device Reporting:

If there is a serious deterioration in the state of your health when using the medical device, report this to your distributor or to the manufacturer and to the competent authorities of your country.

Technical Data

1. EMG

- 1.0 Dual Channel
- 1.1 EMG Range: 0.2 to 2000 μ V RMS (continuous)
- 1.2 Sensitivity: 0.1 μ V RMS
- 1.3 Selectable Bandpass filter - 3db Bandwidth,
 - a. Wide: 18 Hz +/- 4 Hz to 370 Hz +/- 10% - Reading below 235 microvolts
10 Hz +/- 3 Hz to 370 Hz +/- 10% - Reading above 235 microvolts
 - b. Narrow: 100 Hz +/- 5 Hz to 370 Hz +/- 10%
- 1.4 Notch filter: 50 Hz - 33 dbs
- 1.5 Common Mode Rejection Ratio: 130 dbs Minimum @ 50 Hz
- 1.6 Batteries: 4 x AA
- 1.7 Work / Rest periods: 2-99 seconds
- 1.8 Number of Trials: 2-99

2. Neuromuscular Stimulation

- 2.1 Dual Channel
- 2.2 0 - 90 mA into 500 ohm load, 70 mA maximum into 1000 Ohm load, 46 mA maximum into 1500 Ohm load. (Indication only, actual intensity varies depending on condition of skin electrodes).
- 2.3 Type: Constant current
- 2.4 Waveform: Symmetrical, rectangular, bi-phasic with net zero DC current
- 2.5 Pulse width selection: 50-450 μ S
- 2.6 Pulse Rate Selection: 1-100 Hz
- 2.7 Work / Rest periods: 0,5-99 seconds
- 2.8 Time 1 - 99 minutes
- 2.9 Ramp up time: 0.1 - 9.9 seconds
- 2.10 Preset and user programmable treatment programs
- 2.11 Automatic output shut off with detection of open electrode (> 0,5 mA)

The graphics of the waveform and of the intensity depending on resistance can be found on page 92.

Dimensions: L 141,6 mm / W 74 mm / H 40,4 mm - Weight of the device: 275 g (incl. batteries)
Firmware version 3.1B or later

Wireless Link:

RF Frequency Range: 2.4000 GHz to 2.4835 GHz / RF Default Output Power: <= 1,0 mW

RF Maximum Output Power: <= 2,52 mW / Antenna Gain: 0 dBi

Max. environmental conditions for transportation: -25 to +70 degree centigrade, 0-93% atmospheric humidity. Recommended environmental conditions for storage: -25 to +40 degree centigrade, 0-93% atmospheric humidity. The device shall not be stored in direct sunlight. Allow up to one hour for the device to return to normal operating temperature if stored in conditions below +5°C or above +40°C.

Environmental conditions for use: +5 to +40 degree centigrade, 0-93% atmospheric humidity, air pressure 700-1060 hPa. During use the surface temperature of the device can be 41°C or higher, but never exceeds 48°C.

The **DuoBravo** is protected against debris larger than 12.5 mm in diameter from entering the cabinet and is not water resistant and rated IP20.

Explanation of Symbols

	Attention Warning		Manufacturer
	Type BF applied parts		Date of manufacture
	Serial-Number of the product		CE-mark: conforms to essential requirements of the MDD 93/42/EEC
	Keep dry		Not suited for persons with cardiac pacemaker
	Batch code		Read the instructions
	Medical Device		Unique Device Identifier

Safety Adjustments

1. The +buttons become inactive after one minute running time of a stimulation program to avoid an unintentional increase of intensity . This lock is switched off again, by pressing the +button three-times repeatedly. Now an increasing of intensity is possible again for another minute.
2. If the electrical circuit is interrupted during operation ($> 2 \text{ mA}$), for example by losing an electrode or by a broken lead wire, the device stops stimulating immediately.
3. When a phase ends and the next phase contains stimulation, the intensity is automatically set to 8 mA and must be set to the desired setting by the user.
4. The **DuoBravo** features a stimulation circuit failure detection system incorporated into the device and its software. Four pulses with a pulse width and voltage not exceeding $1 \mu\text{S}$ and 6 volts are triggered before each main pulse to check the integrity of the circuit. These pulses cannot be felt and do not effect the efficacy of the device. If a failure is detected during this test, STIM F will display on the LCD. If STIM F is repeatedly displayed on the LCD then the unit should be returned to the manufacturer for testing. This is an added safety protection system to ensure each pulse generated by the device is safe for the user.

Introduction – What is DuoBravo ?

The **DuoBravo** is a modern device for biofeedback and neuromuscular stimulation. The term biofeedback refers to a method to perceive biological processes with the help of a technical device. This instrument measures physical activity - even if there are no visible or no voluntary muscle contractions - and converts them to optical and/or sound signals. **DuoBravo** is a combined two-channel device that can be flexibly used for electromyography (EMG) or EMG-triggered stimulation (ETS) or for neuromuscular stimulation. Doctors and patients are the intended operators for this unit.

The **DuoBravo** is available in two different versions. **DuoBravo U** can be used for urologic, gynecologic and proctologic therapy while **DuoBravo N** is designed for neurologic and orthopedic treatments. Each version includes preset EMG- and ETS-programs as well as multiple stimulation programs in different modes. Additionally the devices have 4 customizable programs in which individual STIM-, EMG- and/or ETS-phases and a combination can be set.

The clip on the back of the unit allows the user to fasten the device to the clothes, the unit is intended to be supported by hand during normal use. The belt clip is for temporary relief only. During treatment it's not necessary to hold **DuoBravo**, since it can stand upright on a flat surface in its stand. This allows the patient to train or to be treated while observing the success of the individual efforts. The muscle activity is indicated by a yellow-green bargraph on each side of the screen and by numerical values in the bottom-line of the display.

At the end of each session, the screen automatically displays the average values of muscle activity in this session. The **DuoBravo** can be used with or without connection to a personal computer. The stored treatment data can be transmitted to the PC at a later date using a separately available documentation tool.

This device is quality controlled and fulfills all legal requirements for medical products. Therefore this device carries the CE-sign „CE 0123“. The serial number is on the rear of the unit.

Manufacturer: MTR+ Vertriebs GmbH
Kamenzer Damm 78, 12249 Berlin, Germany
www.mtrplus.com

Indications

With advice of a physician or physiotherapist, DuoBravo can be used by the patient without any help to support muscle training for:

Main indication for DuoBravo N:

- Stroke Rehabilitation
- Paresis
- EMG-Biofeedback Training
- Relaxation

Side indication for DuoBravo N:

- Improvement of blood flow
- Improvement of muscle spasticity
- Muscle Strengthening
 - Avoidance of muscular imbalances
 - Activation of muscular fibers
 - Muscle Fatigue
 - Strengthening of atrophied muscles
 - Muscle Recovery
- Conservation and extension of mobility
- Release of bodily endorphins
- General Pain Relief and Pain Therapy including but not limited to:
 - Arthritis
 - Back Pain including but not limited to:
 - Chronic Back Pain
 - Cervical Spine Pain
 - Thoracic spine Pain
 - Lumbar spine Pain
 - Circulation problems
 - Fracture Pain
 - Myalgia
 - Neuralgia
 - Phantom Pain

Contraindications

DuoBravo may not be used in the following cases:

- By patients with electronic implants like pumps or a pace maker (as functional disturbances can occur).
- By patients with myocardial damage (damage to the cardiac muscle) or arrhythmia (inconsistent heart beat).
- In presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide.
- When operating machinery with a potential of danger.
- While driving a vehicle.
- By patients with epilepsy.
- If severe local inflammation or thrombosis is present.
- By patients with mechanical ileus.
- By patients with febrile illnesses or infectious diseases.
- Do not use with a surgical high frequency apparatus because burn marks could occur under neath the electrodes.
- When using metal-containing contraceptive spirals.

Please pay unconditional attention to the following basic rules:

- Do not stimulate transcerebraly (via the head).
- Do not stimulate via the Carotis-Sinus-Nerve (carotid artery).
- Do not stimulate in the area of the Sinus knot (heart).
- Do not stimulate on the Bulbous (eye).
- Do not stimulate in the Larynx-Pharynx-Area (throat / mouth).
- Do not stimulate in places with skin defects.
- Do not position the electrodes directly on the spine.

Take special care with the use of the electric stimulation:

- if stimulation is located close to osteosynthesis material (metal in the body).
- if after an operation the healing process can be interrupted through a muscle contraction.
- after acute bleeding or bone fractures with a tendency to bleed.
- during menstruation or pregnancy.
- in case of de-sensitization of the skin (damaged nerves).
- if patients cannot use the device properly, due to mental or physical limitations.

Warnings

- **DuoBravo** is type BF equipment and may not immersed in liquids.
- If using rechargeable batteries be sure to use a CE approved battery charger. Never connect the **DuoBravo** directly to the mains, to a battery charger or to any other mains powered equipment.
- The operator may not touch the battery contacts & patient simultaneously.
- To avoid the effects of electromagnetic interference never use the **DuoBravo** in the EMG-mode, with less than 4 metres distance to a cellular telephone or near any other powerful radio interference producing equipment that causes electrical signals etc. In the EMG Mode, the **DuoBravo** may be susceptible to strong interfering radio type emissions that may lead to temporally increased EMG microvolt readings. If suspected interference is causing a problem with EMG readings, change position of the device and move away from any electrical emitting equipment. The reading will immediately return to the correct value when the interference ceases.
- The power supply of some laptop computers or devices which do not have an earth connection may cause significant interference affecting the DuoBravo. As a precautionary measure please make sure that the power cord of the laptop or the devices are as far away as possible from the lead wires of the patient. If interference persists, remove the power cable from your computer and use it with its battery. Winding the lead wires around a limb (arm or leg) may also help this interference.
- If device operation occurs close to short wave, microwave, cellular or other wireless equipment, instability of the values of the stimulator could occur causing a temporary reduction in current. Move away from the source of interference and continue the program.
- Burn marks could occur underneath the electrodes if used with surgical high frequency apparatus.
- The application of the electrodes near the thorax can increase the risk of cardiac fibrillation.
- Patients with metal or electronic implants may only stimulate with permission of an experienced doctor.
- Stimulation in the area of tumors is allowed only by instruction of an experienced doctor.
- If high frequency is combined with a wide pulsewidth, please stimulate with low intensity only.
- If you have swallowed any part of the device or any of its accessories, please contact a doctor immediately.
- TENS and muscle stimulation should be used during menstruation or by pregnant women as well as by children under the guidance of a doctor or therapist only.
- The **DuoBravo** may not be operated if the battery compartment is not closed. Before opening the battery compartment, make sure that there are no lead wires in the jacks at the head of the unit.
- Remove batteries and battery cover, drain water from device and allow (DuoBravo) to dry normally for 48 hours if accidentally submerged in water.

Please, keep this device stored out of the reach of children. Do not leave the lead wires in reach of children, since strangulation is possible due to the length of the lead wires.

Advice on Application of Electrodes

Make sure that **DuoBravo** is turned off (no indication visible on the LCD display).

Now connect the single poled ends of the lead wires with the electrodes.

The output plugs of the electrode lead wires will be plugged into the channels A or B of the **DuoBravo**. The reference lead wire will be connected to its jack as well if EMG will be used. Furthermore make sure, that all plugs fit completely into the electrode sockets as well as the device sockets.

Please don't use oil, ointment or cream, since all these have an isolating effect.

Now place the electrodes according to the instructions of your doctor or physiotherapist. The self-adhesive skin electrodes are provided with a sticky gel and can be re-used after removal from the skin. Replace these electrodes when the adhesive force is weak.

Electrodes which are damaged due to careless handling may not be used.

Electrodes with damaged lead wire insulation (e.g. visible copper wire) may not be used.

Electrodes with damaged, missing or stained aluminium coating or carbon layer, also with insufficient adhesive gel may not be used.

Electrodes may not be altered or tampered with (e.g. by scissors or other tools).

Do not pull the electrodes by force on the electrode lead wire. To remove the electrodes from the skin, please lift the electrode at the edge towards the middle.

Do not use electrodes with ointments or creams (especially which enhance the blood flow or act as local anaesthesia).

Apply the electrodes on clean and healthy skin only.

In the case of strong bodily hair growth, the hair should be shaved or shortened with scissors.

Neither stretch the skin nor the electrodes when applying the electrodes.

Effective current densities of more than 2 mA/cm² require an increased attention of the user with regard to the intensity adjustment.

For standard treatment we recommend the use of our self-adhesive electrodes with a minimum area of 40x40 mm. This size is suited for the reference electrode as well. For best results please make sure that it will be placed on a bone near the treatment area, not on muscles.

After all electrodes have been removed, wash the skin with mild soap and rinse thoroughly. Then dry your skin.

All electrodes are for single person use only.

Skin electrodes can contain amounts of nickel. Allergy sufferers should check this before using.

Should skin irritation occur whilst using the **DuoBravo**, please interrupt the application and inform your doctor, so that an allergic reaction can be excluded.

Maintenance

DuoBravo has been designed to be maintenance free. Nevertheless, please adhere to the following advices:

If needed, please clean **DuoBravo** with a mild detergent, water and a damp cloth. Do not expose **DuoBravo** to extreme humidity and do not submerge under water. You may clean the lead wires with a damp cloth. You may also from time to time treat the lead wires with talcum powder in order to prevent cracks and to improve the lead wires life. A clean cloth, moistened with an alcohol free disinfection liquid can be used to clean the device. Please follow all instructions of the manufacturer of the liquid.

After termination of the therapy, **DuoBravo** can be re-used by other patients. Please clean the device with a mild disinfectant cleaner.

The **DuoBravo** does not require calibration. The **DuoBravo** is not repairable and contains no user serveable parts. No modification to the **DuoBravo** is allowed.

Warranty

MTR+ Vertriebs GmbH provides a 24 month warranty, commencing from the date of our delivery (invoice date), that in case of failure, the **DuoBravo** will be serviced or replaced free of charge. All returns must first be authorised by MTR+ Vertriebs GmbH in advance.

Excluded from this warranty are failures of the **DuoBravo** which result from excessive use, careless or improper handling as well as damage by force.


The warranty loses its lawful meaning if the producer MTR+ Vertriebs GmbH has not given its consent in writing prior to repairworks or servicing being carried out. Without this written consent the manufacturer will not return the device.

The manufacturer reserves the right of technical changes, which may serve to increase treatment success, the functionality of the device or its accessories.

Operating Elements

The **START**-button allows the choice of automatic (press) or manual (hold the button for 3 seconds) threshold-setting if an EMG-or ETS-program is chosen. The program will start after the **START**-button is released.

+/- buttons for threshold-setting and setting of intensity for channel A, as well as for selection of parameters in the programming-mode

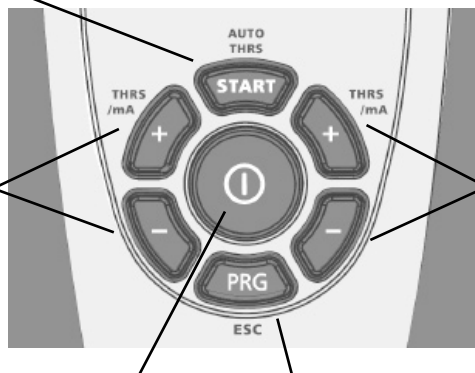
+/- buttons for threshold-setting and setting of intensity for channel B, as well as for altering the parameters in the programming-mode

Central On-/Off-button

PRG-button for choice of program and for storing the individual parameters. Hold PRG for one second to select next phase.

Battery Compartment


LOCK-button for fixing one single program and for uploading statistics to a PC

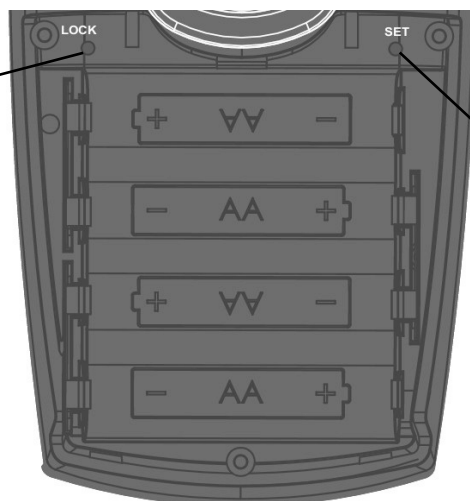SET-button for basic settings and configuration of custom EMG and ETS programs

Scheme for placement of 4 batteries type AA

Connections

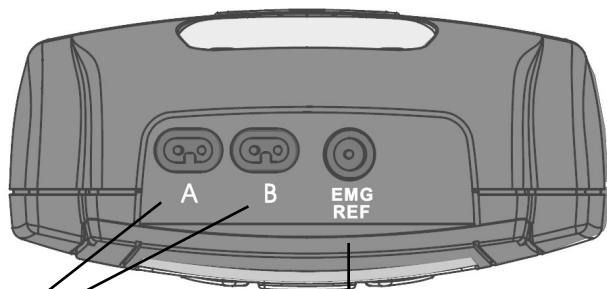
Jacks to connect the electrode-lead wires of channel A and/or B, each suited for surface-electrodes.

Jack for the single pin EMG-reference lead wire (skin-electrode)

The Display

After selecting a program with the **PRG**-button **DuoBravo** shows information about those parameters in its “Home Screen” that are used in the first phase of the program (stimulation frequency, pulse width, number of trials, threshold, etc.). Start the treatment by pressing or holding **START**, the display changes to work-screen and shows necessary information for treatment in the active phase.

Work-Screen (in EMG-mode)

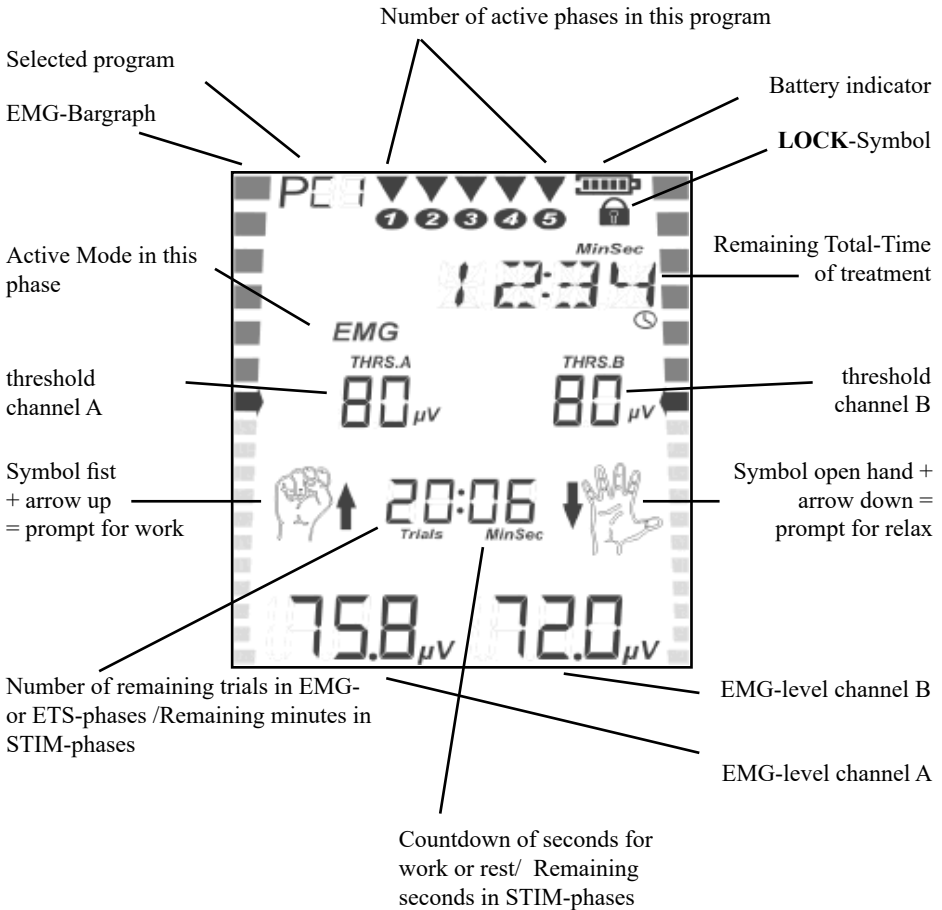The Info-Line

The info line in the upper part of the screen shows always the remaining time of the entire program during treatment.

Outside a treatment session, it provides more additional information:

- When selecting a program, the target of the treatment is indicated in abbreviations. (REHAB = Rehabilitation, etc.).
- If the circuit is interrupted during stimulation, caused by removal of an electrode or a broken lead wire, the term CABLE flashes in the info line.
- It shows custom values to be set.

The info line displays all custom settings for configuring the device. In the following you find all expressions and abbreviations that are used in the programming-loop in alphabetical order:

CH.B. ON / CH.B. OFF: Channel B is switched on or off

DELAY: Possibly delay of current between channel A and B if synchronous stimulation is selected

EMG OFF / EMG ON: EMG-measuring during rest-period of STIM-phases

FB. CONT: Continuous acoustic feedback

FB. OFF: No acoustic feedback

FB.ABOVE / FB.BELOW: Feedback above or below the threshold

FREQ.: Number of stimulation pulses per second

INHIB.OF: Inhibition switched off

INHIB.ON: Inhibition of feedback on channel B

NARROW: Narrow bandpass filter used

PC OFF: Connection to PC is switched off

PC SWC: Connection to PC by secured wireless connection

PULSE: Duration of each single pulse in μ sec

RAMP DN.: Time needed from the end of complete intensity to complete end of current

RAMP UP: Time to reach complete intensity of stimulation

REST: Duration of rest-sequences in seconds

STIM.TME.: Duration of EMG triggered stimulation in seconds

THRESH: Threshold-adjustment during treatment automatic (AUTO) or manually (MAN)

TIME: Duration of a STIM-phase in minutes

TRIALS: Number of repetitions

TRIGGER: Time that the threshold must be reached to trigger stimulation

VOLUME: Sound volume of beeps and spoken instructions

WIDE: Wide bandpass filter used

WORK: Duration of work-sequences in seconds

Basic Settings

In all programs, even in the preset programs, individual settings can be made as follows. After you have put in four batteries type AA (s.p. 65) and you have turned on the device with the **ON/OFF**-button, you enter the configuration mode by pressing the **SET** button in the battery compartment for three seconds. Now you can adjust the settings for the first phase of the selected program. To configure the next phase, **SET** should be pressed again.

The parameter that can be changed is shown in the info line and can be altered by pressing the + button of channel B. Pressing the + button of channel A will confirm the choice and the next variable parameter is displayed. This system - changing settings with B / storage and stepping forward with A - remains always the same during the entire process.

Please note that the following settings you make here, are automatically stored for all subsequent phases and all other programs as well. Of course you can change any setting at any time in the same way again.

STIM-phases

- Channel B can be switched on or off.
- Work Rest phases have these options:
 - a) Turn on / off measurement of EMG during rest period.
 - b) Selection of wide or narrow filter.
- PC - Select off or SWC (wireless).

EMG-phases or ETS-phases

- Sound-volume
- Automatic threshold-adjustment on or off
- Wide- or narrow-filter for channel A
- Channel B switched on or off
- Filter for channel B if it is in use
- PC - Select off or SWC (wireless)

In the customizable programs PC1 - PC4 of course you have the ability to set additional parameters in the same way as described on pages 66-69.

Setup Real Time Clock

Before first use of **DuoBravo** the integrated “Real Time Clock” should be set. It’s necessary later for therapy-control. **DuoBravo** has a small battery that will support the clock with energy for several years. For this reason, there is no need to set again, even if the 4 main batteries have been substituted or removed. But we recommend to check and re-adjust periodically.

To set time and date, please select any program and press the **SET**-button for at least 3 seconds. You’ll reach the configuration-mode. Now press **SET** once again for 10 seconds. On screen HOUR is indicated and you can set the hour with the +button of channel B (00-23).

By pressing the +button of channel A you confirm this choice and step forward to the setting of MINUTE. This is indicated on the screen as well, giving the opportunity to select the minutes via +button of channel B in the same way as before.

This way you are able to set one after the other:

- HOUR
- MINUTE
- DAY
- MONTH and
- YEAR

If all settings are finished, please press the **ESC**-button on the frontside of the device. All values will be automatically stored and the clock starts working.

Setup Language

DuoBravo is able to show 6 different languages. If you are in the setting-mode of the RTC (see above), please press the **SET**-button for 10 new seconds once again. Now you can select German-, English-, French-, Spanish-, Italian- or Polish-language with the B-buttons. Finally you store your choice with the **ESC**-button.

Preparation of Therapy

1. Open the battery compartment at the back by pushing down the lid and insert 4 batteries type AA by hand according to the instructions on page 51. Then shut the battery-compartment.
2. Please make sure, that **DuoBravo** is switched off (no indication in LCD-display). Connect the single pin ends of the electrodes plugs with the surface-electrodes.
3. Insert one electrode-lead wire into the jack named A in the head of the device. If you want to use both channels, insert the second lead wire into jack B.
4. Additionally, please put the black round EMG-reference plug into the bigger round jack, named "EMG-REF" and connect the other ends to a skin-electrode.
5. Furthermore, make sure that all plugs fit completely into the electrode sockets as well as the device sockets.
6. Now place all electrodes according to the instructions of your therapist on clean and healthy skin only.
7. Now switch on **DuoBravo** by pressing the central **ON/OFF**-button once. You can see the version of the device for 3 seconds before the last used program appears, showing its parameters.
8. If the device is not locked to a single program, you can select your kind of treatment by repeatedly pressing the **PRG**-button. A great variety of programs with up to five phases each are available. Within one single program the phases can be programmed in different modes.
9. After selecting a program, you always see the "Home Screen" first, with all information regarding the parameters, that are contained in the first phase of this program. If the chosen program has an EMG- or ETS-phase, the manual threshold-setting is possible. In this case, the treatment starts by holding the **START**-button for more than 1 second. If you prefer the automatic threshold-setting, no adjustments have to be done and treatment starts by short pressure of the **START**-button. Follow the work/rest prompts, so the device can automatically select the threshold setting for you.
10. EMG is susceptible for interference. For best results do not use lead wires longer than 1.5 m and keep the lead wires as close to your body as possible.

The Different Phases

DuoBravo offers programs that can consist of one or more phases. There are basically three different types of these phases which can be therapeutically combined in a useful way.

a) **EMG** (= Electromyography) - purely active biofeedback training is carried out. The patient tightens the muscle according to the requirements in the display for the specified time and then releases it again. The device gives an acoustic and visual feedback about the strength and duration of contraction and relaxation. If the relaxation is target of the treatment, the positive signal will be given in the active/work period (that is indicated by the fist) as well, if muscle tension will be below the threshold. In an EMG-phase the device only measures muscle activity and doesn't supply any current.

b) **STIM** - In this phase the muscle is stimulated passively.

- Continuous current (CON), which does not include rest periods. The parameters do not alter.
- Modulated current, which also flows without a rest, but frequency and pulse width change in a predetermined rhythm. This creates a feeling of rising and falling intensity.
- Current in work/rest-rhythm, where a muscle is stimulated in each work period to cause contraction. In all subsequent rest periods the muscle shall be relaxed as far as possible, to avoid strain.

The effect of current delivered by the device will vary depending on configuration. It may have circulatory-enhancing, detoxifying, relaxing, strengthening and pain-reducing effect. During the stimulation **DuoBravo** cannot measure EMG, but the muscle tone can be measured in the relaxation sequences (the rest-period). This control may be necessary to avoid any overstrain of the treated muscle.

c) **ETS** (= EMG-triggered stimulation) - this is a combination of active biofeedback training (see above - EMG) and immediately afterwards triggering passive electrical stimulation (see above - STIM). If the patients succeed to reach a preset threshold for the muscle tension on their own, the device automatically transmits current that causes the muscle to contract. The movement generated by the current, is initiated by the brain using this method. It's the target to improve the force of one's own muscle contraction. In neurological treatment, it is often intended to remove interferences in neural pathways, so that the muscle can be activated without the help of a device.

Therapy with DuoBravo

DuoBravo is a versatile device, that measures, compares, displays results, adapts itself, gives current, calculates statistics ... and much more.

For this reason, each kind of therapy is possible without complicated preparation. Besides four customizable programs, the following number of pre-set programs are available, that can be used for therapy at once, without any further adjustment.

- **DuoBravo N** (Neurology) = 12 programs

Use of pre-set programs

1. Connect the electrodes to the device (s.p. 83)
2. Place electrodes according to the instructions of your therapist
3. Switch on **DuoBravo** with the central **ON/OFF**-button
4. Select the desired program with the **PRG**-button
5. Press **START**
 - In phases with passive stimulation you'll hear a short melody and the mA segments will flash, showing that the current can be adjusted with the +button of any used channel. The intensity will increase 1 mA with each press - please adjust current to the desired (comfortable) level ensuring no pain is felt. Please note that this increase is only possible, if the +button is active (s.p. 58 - Safety Adjustments).
 - In EMG- or ETS-phases threshold is set automatically by performing two short work/rest-periods. If that's not required you can switch off this function by holding the **START**-button for three seconds. The threshold can be manually set with the A+ and B+ buttons.
 - An EMG-phase will start with the first rest-period now. If you are using an ETS-phase the short melody sounds and flashing mA prompt for intensity-adjustment.

If the selected program has several phases, each phase-end will be indicated by signal beeps. Please note that stimulation-intensity will fall down to 8 mA each time a new phase starts. The intensity will need to be adjusted accordingly.

Therapy with DuoBravo

Therapy in a phase with stimulation

In phases with work/rest-rhythm you can support treatment success by active muscle contraction during work-period effectively. Please take care that the muscle will relax in rest-periods completely. In display **DuoBravo** shows the fist as symbol for “work” and the open hand for “rest”.

At the beginning of each stimulation phase you hear a short melody and the mA-indicator is flashing in the display, showing that intensity should be adjusted now with the +/- buttons of each active channel.

The stimulation intensity is indicated with numbers in the bottom-line of the display and is shown by the lateral bar graphs on the screen as well. The more bars are active, the higher the current intensity.

If you want to interrupt the treatment of a stimulation-phase, please press **START** and the current will stop immediately. In the infoline of the display the term PAUSE appears. You can continue the treatment by pressing **START** once more. Then stimulation starts again for the remaining treatment-time with the same intensity it had before it was paused.

In the rest-periods of programs with work/rest-rhythm, EMG can be measured to see if the treated muscle is relaxed enough to avoid overloading (s.p. 69 and 78). The corresponding values are shown in the bottom line of the display.

Therapy with DuoBravo

Therapy in an EMG-phase

The device indicates a pre-set rhythm of work and rest, that the patient shall execute as precisely as possible.

During work-period the symbol of the fist with up arrow is seen on the screen. The feedback regarding the success of your own muscular activities is given by the increasing numbers in the lowest line and by the yellow/green bargraph on each side of the display. Reaching the green area shows that the threshold - regardless whether it has been set automatically or has been set manually - is met and the trial was successful. Additionally this success is indicated by beeps. Please take care to keep your efforts active, until the time, that is indicated beside the fist, has counted down to zero and the symbol of the open hand with arrow down appears. Then relax the muscle as much as possible.

If you have chosen automatic threshold adjustment during treatment, threshold would be reduced if you are not successful and it would be easier to succeed in the next trial. If this next trial is successful again, the threshold will increase.

The process of work and rest shall be repeated until the number of remaining trials counts down to 00 (s. p. 67). Please note that the display doesn't change before the relaxation-phase is completed.

If relaxation is target of the treatment, the positive signal will be given during the work period (fist) as well, if muscle tension is below the threshold.

Therapy with DuoBravo

Therapy in an ETS-phase

An ETS-phase is similar to an EMG-phase, which means, that you follow the rhythm of work and rest, given by the device, with your muscular activities. If you succeed reaching the threshold the term BRAVO is indicated in the infoline of the display and the device stimulates your muscles automatically. In this case your muscle will be contracted passively. Please support the treatment success by keeping your muscle contracted during stimulation until there is no current-flow anymore. Then the rest-period starts and you should relax the muscle as much as possible.

If you don't reach the threshold in one of the trials, you will not receive stimulation and the rest-period starts at once after termination of the work-period.

If you have chosen automatic threshold adjustment during treatment (s.p. 83), threshold will be reduced in this case and it would be easier to succeed in the next trial. If this next trial is successful, the threshold will increase.

The process of work and rest shall be repeated until the number of remaining trials counts down to 00 (s. p. 67). Please note that the display doesn't change until the relaxation phase is completed.

Therapy in combined programs

Programs, which contain several modes (STIM and/or EMG and/or ETS), will always start with automatic threshold-setting, unless you skip it (s.p. 82). When threshold is set you can adjust the current-intensity with the +/--buttons for channels A and(or) B if the program contains a STIM-phase. After this, each step is identical to the description above.

Additional Options for Configuration of DuoBravo

Additionally to the pre-set programs as described in the previous chapter, **DuoBravo** offers a great variety of possibilities to adapt the therapy-parameters according to individual needs.

This concerns the pre-set programs that can vary in regard of the basic adjustments as well as the 4 customizable programs, in which the experienced therapist or user is able to create completely new programs with up to 5 different phases.

Configuration of Custom Programs

Creating a new program is not difficult, because the programming always follows the same easy scheme. The device leads in a loop to each parameter that can be chosen - indicated by flashing bars.

1. Switch on **DuoBravo** with the **ON/OFF**-button
2. Select PC1, PC2, PC3 or PC4 with the **PRG**-button. In **DuoBravo N** these custom programs are 13-16. In the left upper corner of the screen PC1-PC4 will be indicated.
3. Keep the **SET**-button in the battery-compartment pressed for 3-5 seconds. Then **DuoBravo** changes to programming mode and the bars beside the first parameter that can be set are flashing.
4. In the following the right “+/-”buttons (channel B) are responsible for changing the flashing parameter, while the left “+” button (channel A) is used to store the choice and step forward to the next parameter. With the left “-”button it is possible to go one step back.

With the described scheme “right side changing / left side storing and step forward” all necessary decisions can be met. If the first parameter is shown a second time, the configuration of this phase is complete and by short pressing of **SET** the next phase is available. If you store “OFF” then, this and all following phases are deleted. If all phases are ready programmed, the program is stored by pressing the **ESC** button. This configuration will stay available until a new program is entered and saved or if the device is reset to its original settings

Parameters of Custom Programs

Programming PC1

PC1 is a multi-phase program for passive stimulation only. With A+ the following parameters can be chosen and then altered by B+, if necessary.

- **Type of stimulation** - (indicated in display: CON / SYN / ALT) - sets whether stimulation shall be continuously or in work/rest-rhythm.

If **CON** (= continuous current) is selected, you can afterwards only decide about:

- **Duration of phase** (PH.TIME)
- **Frequency** (FREQU.) - fixes the number of pulses or pulse-groups per second. Frequency is indicated with Hz.
- **Pulse width** (PULSE) - names the duration of each single pulse and is measured in microseconds (μ S).
- **Channel B** (CH.B. ON / CH.B. OFF) switched on or off
- and **PC-Interface** (PC SWC / PC OFF) - here you can decide whether all treatment data will directly be transferred during treatment to a PC with optional documentation tool by secured wireless connection or not. During home-therapy PC OFF should always be selected to save battery-power.

If you want to use both channels in a work/rest-rhythm simultaneously, please select **SYN**. In this case the following parameters can be set:

- **Duration of phase** (indicated in display s.a.)
- **Frequency** (s.a.)
- **Pulse width** (s.a.)
- **Duration of work-period** (WORK) - sets the time that one single muscle contraction is held for
- **Duration of rest-period** (REST) - sets the time of relaxation after contraction
- **Ramp up-time** (RAMP UP) - means the time (in secs.) of increasing current
- **Ramp down-time** (RAMP DN) - means the time (in secs.) of decreasing current
- **Channel B** (CH.B. ON / CH.B. OFF) switched on or off
- **Delay of channel B** (DELAY) - only if needed
- **Measuring of EMG in rest-periods** (EMG OFF / WIDE / NARROW)
- and the **PC-Interface** (s.a.)

ALT should be selected if alternate muscle contractions are required (i.e. channel A contracts followed by Channel B). All parameters are set identically to SYN mode excluding delay which is not available in this mode.

If programming of one phase is complete, press **SET** to go to the next phase.

The range of settings and the size of possible steps is shown on page 69.

Parameters of Custom Programs

Programming PC2

PC2 is a single-phase EMG-program. With A+ the following parameters can be chosen and then altered by B+, if necessary.

- **Volume** - (indicated in display: VOLUME) - sets the volume of sounds for feedback in steps 1-10
- **Duration of work-period (WORK)** - sets the time one single muscle contraction is held for
- **Duration of rest-period (REST)** - sets the time of relaxation after contraction
- **Number of trials (TRIALS)** - sets the number of active exercises that shall be done
- **Threshold (THRESH)** - defines whether automatic threshold-adjustment during therapy is switched on (AUTO) or off (MAN).
- **EMG Sound Feedback** – a) Off (FB.OFF). On (FB.ON), Above Threshold (FB.ABOVE) – normally used when training muscle force, Below Threshold (FB.ABOVE) – normally used for relaxation training.
- **Wide or narrow filter (WIDE / NARROW)** for channel A - this choice is necessary to avoid wrong measurement, caused by heart beat (60-120 Hz). Narrow filter should be used for the torso, wide filter should be used for legs and arms.
- **Channel B (CH.B. ON / CH.B. OFF)** switched on or off
- **Wide or narrow filter** for channel B
- **Feedback inhibition (INHIB.OF / INHIB.ON)** - will be selected, if the user shall get positive feedback only, if (s)he will be above threshold of channel A **and** below threshold of channel B. This method is used to avoid wrong EMG-results, caused by contraction of wrong muscles, which might be near the muscles that shall be measured.
- and **PC-Interface (PC SWC / PC OFF)** - here you can decide whether all treatment data will directly be transferred during treatment to a PC with optional documentation tool by secured wireless connection or not. During home-therapy PC OFF should always be selected to save battery-power.

Parameters of Custom Programs

Programming PC3

PC3 is a single-phase ETS-program. With A+ the following parameters can be chosen and then altered by B+, if necessary.

- **Volume** - (indicated in display: VOLUME) - sets the volume of sounds for feedback in steps 1-10.
- **Duration of work-period (WORK)** - sets the time that one single muscle contraction is held for.
- **Duration of rest-period (REST)** - sets the time of relaxation after contraction.
- **Number of trials (TRIALS)** - sets the number of active exercises that shall be done.
- **Duration of stimulation-period (STIM.TME)** - sets the time of current-flow if threshold is reached.
- **Frequency (FREQU.)** - sets the number of pulses or pulse-groups per second. Frequency is indicated with Hz.
- **Pulse width (PULSE)** - sets the duration of each single pulse and is measured in microseconds (μ S).
- **Ramp up-time (RAMP UP)** - sets the current rising time.
- **Ramp down-time (RAMP DN)** - sets the current decreasing time.
- **Threshold (THRESH)** - defines whether automatic threshold-adjustment during therapy is switched on (AUTO) or off (MAN).
- **EMG Sound Feedback** – a) Off (FB.OFF). On (FB.ON), Above Threshold (FB.ABOVE) – normally used when training muscle force, Below Threshold (FB.ABOVE) – normally used for relaxation training.
- **Wide or narrow filter (WIDE / NARROW)** for channel A - this choice is necessary to avoid wrong measurement, caused by heart beat. If the patient measures EMG near the heart you should use NARROW-filter.
- **Channel B (CH.B. ON / CH.B. OFF)** switched on or off
- **Wide or narrow filter** for channel B
- **Trigger-time (TRIGGER)** - is the time that muscle-tension of a user must remain about the threshold in order to trigger the current. There are two settings available. 0.25 seconds is suitable for urology and avoid excess muscle twitches triggering the current. The patient must hold the EMG value above the threshold for 0.25 seconds or more to trigger the current. 0.01 seconds is suitable for neurological treatments and can detect the smallest muscle twitch which is above the threshold. The patient only needs to hold the EMG value above the threshold for 0.01 seconds to trigger the current.
- and **PC-Interface (PC SWC / PC OFF)** - here you can decide whether all treatment data will directly be transferred during treatment to a PC with optional documentation tool by secured wireless connection or not. During home-therapy PC OFF should always be selected to save battery-power.

Parameters of Custom Programs

Programming PC4

PC4 is a multi-phase program, in which each of the max. 5 phases can be configured in STIM-, EMG- or ETS-mode. Press SET to select consecutive phases.

If you want to change a configured program and you do not need all the phases, please set the first phase that is no longer needed to "OFF". This and all further phases are deleted.

Please note that for safety reasons each phase containing stimulation will start with 8 mA max. Intensity has to be adjusted newly to a therapeutic effective level without causing any pain.

The settings can be made as follows.

- Frequency: 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20 Hz, then increasing in 5 Hz steps until 50 Hz and in steps of 10 Hz up to max. 100 Hz.
- Pulse width: 50 to 450 μ S in steps of 10 microseconds
- Work time: 1-99 seconds in steps of 1 second
- Ascent time: 0.1-9.9 seconds in steps of 0.1 second (Ramp up)
- Descent time: 0.1-9.9 seconds in steps of 0.1 second (Ramp down)
- Rest time: 1-99 seconds in steps of 1 second
- Delay: 0.0-4.0 seconds in steps of 0.1 second
- Session time: 1-99 minutes in steps of 1 minute

Lock-Button

If required, the therapist can lock out all of programs.

To do this, select the desired program and press the **LOCK**-button with a thin pin. For this you can use the electrode connector. You'll hear a short beep and the symbol of a padlock appears on screen.

If you want to use any other program later, or if you need to unlock all programs, **LOCK** must be pressed in Home Screen once again and CONN-PC (= PC-connection) or NO.STATS (= no statistics) will appear. Keep the **PRG** button pressed for 3 seconds, the display changes to Home Screen and all programs are available. Press **PRG** for less than one second to return to the Home Screen with the unit still locked.

Threshold

The threshold is a target value that shall be achieved by personal effort. It is measured in microvolts and displayed for both channels in the middle line of the screen. The arrows in the middle of the bargraph on the left and right side of the display represent the thresholds. Since all programs with EMG or ETS phases refer to a threshold, it is necessary to set the threshold for all programs with EMG or ETS phases.

Setting of threshold before start of treatment

By default, the threshold is determined automatically when **START** is pressed if any phase in a program is EMG or ETS.

Pressing the **START** button causes the appearance of the open hand in the display, while the time-indicator for this rest-period runs to 0. At 0, the display switches to the symbol of the fist and the time for this 5-second period of work runs to zero again. **In this short work period, you should try to contract the muscle with the maximum possible effort.** Subsequently a rest period follows and a further work-period. During these two initial work/rest trials for threshold-detection the word TEST is indicated in the infoline. The threshold is calculated automatically after the second phase of work, and is set at 80% of the average value of the period with higher performance. Please note that the first second of each work-phase is not relevant for the average value.

Now the treatment begins with the threshold set and the total time counts down in the infoline.

Alternatively, the automatic setting can be skipped. This is achieved by holding the **START** button pressed for more than one second. In this case, the threshold of the last treatment appears and could be altered manually by pressing the +/- buttons before starting the program.

Threshold

Automatic adjustment of threshold during treatment

Once the threshold is set, it will stay the same for the entire duration of the treatment until you make a new setting - except if you have selected “automatic threshold adjustment.” In this case, the threshold would be adjusted during the treatment according to your performance.


EMG-Phase

In an EMG-Phase the new threshold will be set to 80% of the average work that have been done in the last trial before.

ETS-Phase

If you reach the threshold during a trial, it increases for the next attempt and if you don't succeed to reach the threshold, it goes down slightly. The stepsize of increase or decrease depends on the time that was needed to reach the threshold in the last trial before.

The following graph will illustrate the principle of adaptation in an ETS-phase.

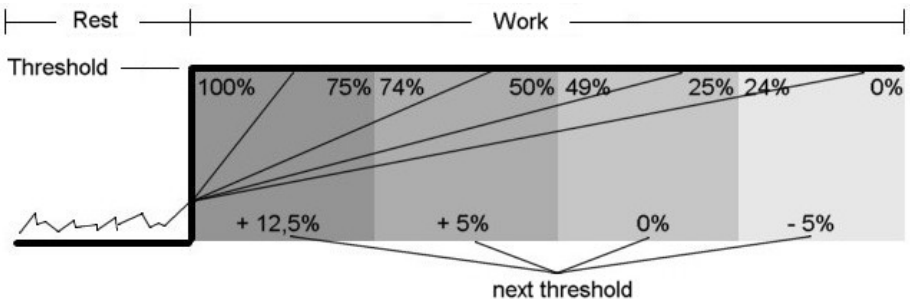If the threshold is reached in the first quarter of the time, it increases by 12.5%.

If the threshold is reached in the second quarter of the time, it increases by 5%.

If the threshold is reached only in the third quarter of the time, it does not increase or decrease.

If the threshold is reached only in the last quarter of the time, it decreases by 5%.

If the threshold is not reached, it decreases by 12.5%.


Therapy-Control

DuoBravo displays statistics after each program has been completed. All values of this completed session, which are named in pages 74-75 will be indicated on the screen automatically. With A+ you can scroll through the statistics. **DuoBravo** saves these statistics to its memory for download to the optional documentation tool.

Up to 200 days of treatment data is stored in the **DuoBravo** memory. When 201 days of data are reached, consecutive days will overwrite the earliest days. For this reason the statistics will always show the actual data of the last 200 treatment days. To access the statistics, press the LOCK button in the battery compartment once if the device was locked in a program - or twice if it was not. The **DuoBravo** displays the number of treatment days saved, the total usage time and the most used program.

Press the A+ button to access the individual treatment days (1 to 200) with the respective number of sessions displayed. Press the START button to check the session data for the selected day. The duration of the session and start time will be displayed. Press B+ to view the start time for up to 10 consecutive sessions for that day. Press the START button again to view the accumulated daily data. Press B+ to scroll through the average daily data for EMG, ETS and STIM. Up to 10 values can be viewed as described on pages 74-75.

Pressing A+ will select the next day

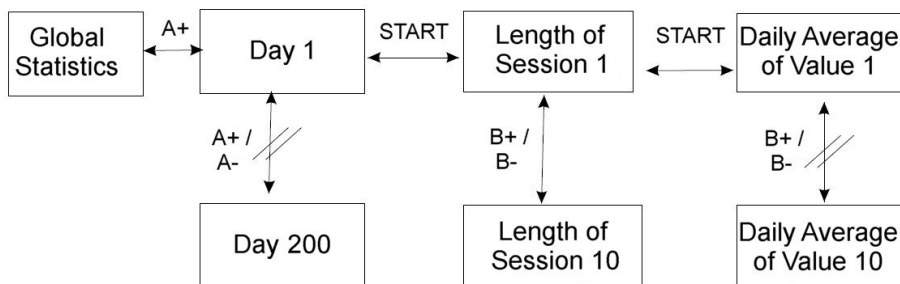Therapy-Control

The additional documentation tool is required to transfer the stored data to a PC. An electronic patient record, can easily be created. This optional documentation tool is very helpful for therapists for statistical reports. To transfer the data please insert the USB dongle, which is supplied with the documentation tool, into the computer.

Please note that data transfer is not possible while checking the daily statistics on the unit.

The transfer of stored data can be done as follows:


- Start the optional documentation tool on your PC.
- Create or select the patient.
- Switch on the **DuoBravo**.
- Press the LOCK-button once, if the device was locked – or twice if it was unlocked.
- The arrows in the display of the unit will flash one after the other indicating the statistics can be transferred to the computer. The statistics will be transferred to the computer when the screen displays the number of days and time used. “NO-STATS” will be displayed if the memory is empty.
- The data transfer starts automatically and your PC shows a message that it is receiving the statistics.
- The transfer may take several minutes depending on the amount of data stored in the **DuoBravo** memory.
- When the statistics have been uploaded, the documentation tool will confirm that the transfer is complete.

The values in the therapeutic memory are retained after they are read, so it is possible to check therapy at any time, without losing information.

If the memory needs to be deleted, it can be done with the following button-combination. Press the LOCK button in the battery compartment until CONN-PC is displayed and then successively SET, B +, B- and SET once again. You will hear a sound to confirm the statistics are being deleted. This action will also reset all custom settings in the customizable programs back to the original state.

Ensure you delete the statistics when new patients use the device.

Statistics

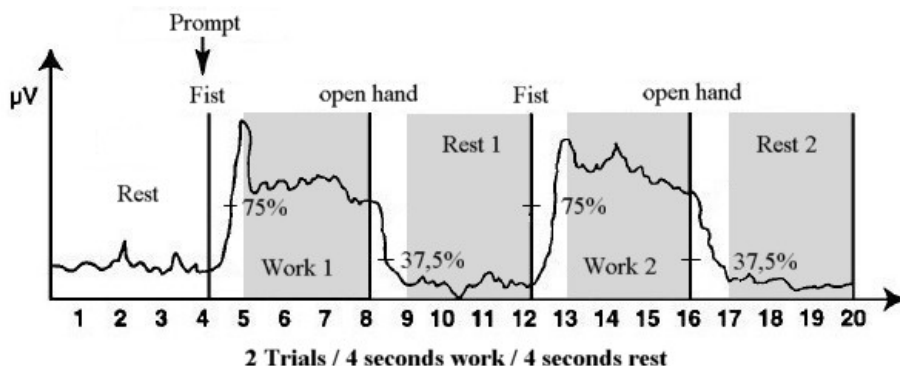After termination of a treatment session **DuoBravo** shows statistics on its screen as follows.

1) WORK.AVG - Work Average:

The average value in microvolts of all the work-segments excluding the first second of each segment.

2) REST.AVG - Rest Average:

The average value of the rest-segments excluding the first second.

3) ONSET - Average Onset:

The average time taken after each respective “work” prompt to reach 75% of the average value of the previous work period. If the onset of trial(s) is longer than 2 seconds, it will be rejected. The display will indicate the average of only those trials, where onset lasts 2 seconds or less. The number of valid onsets is displayed above the onset value. If the number of valid onsets is lower than the number of trials, it means some onsets were rejected.

4) RELEASE - Average Release:

The average time taken, after the “rest” prompt to reach 37.5% of the average EMG from the previous work period. If any release of the trials is longer than 2 seconds then it will be rejected. The display will indicate the average of only those trials where release lasts 2 seconds or less. The number of valid releases is displayed above the release value. If the number of valid releases is lower than number of trials, it means some releases were rejected. If all the releases are above 2 seconds, ABOVE2 will appear on the screen.

5) W/R PEAK - Work/Rest Peak:

It is the peak value during the whole cycle of all the trials.

Statistics

6) WORK DEV - Average deviation (per second) of work:

It is the average deviation of the quarter-second samples from the average value in the second in which they occur (excluding the first second of each work-segment). If deviation is more than 20% it is considered problematic.

7) REST DEV - Average deviation (per second) of rest.

The average deviation of the quarter-second samples from the average value in the second in which they occur (excluding the first second of each rest-segment). If deviation is more than 20% it is considered problematic.

8) ETS AVG - Average intensity of current:

Average intensity of current in all phases with ETS.

9) SCORE - Percentage of success:

The SCORE indicates the speed which the patient can activate the stimulation for ETS. Values closer to 100% indicate a faster average response time and better muscle condition.

10) STIM AVG - Average intensity of current:

Average intensity of current in all phases with pure stimulation.

Trouble Shooting Guide

Fault	Possible Cause	Solution
Trip off - Electrode symbol flashing, current won't increase above 8 mA)	Electrode lost adhesion to skin.	Check electrodes and replace if necessary.
	Electrode too dry	Replace electrode.
	Lead wire not connected properly.	Check lead wire connections and ensure they are fully inserted.
	Lead wire is faulty.	Replace lead wire.
	Device is damaged.	Return device to manufacturer.
Device does not power on.	Flat battery or no batteries installed in device.	Check / Replace batteries.
	Battery installed in wrong orientation.	Remove batteries and install according to the advice on page 53
	Battery connector damaged or flattened.	Remove batteries, check metal parts in battery compartment for damage.
	Device is damaged.	Return device to manufacturer.
Device powers off during use.	Batteries almost flat.	Check battery icon during use and replace batteries if low.
	Device is damaged.	Return device to manufacturer.
Current output is unstable	Device is in use near a device emitting strong radio or magnetic interference.	Move at least four meters away from the offending device and resume treatment.
EMG reading is unstable	Reference electrode is missing	Place reference electrode.
	Device is being used with a notebook computer.	Disconnect power from notebook computer if being used.
	Device is in use near a device emitting strong radio or magnetic interference.	Move at least four meters away from offending device.
	Output peaks in conjunction with heart beat while electrodes are attached near the heart.	Device is reading the heart beat. Change device to Narrow band to filter out heart beat.
Treatment is uncomfortable when using skin electrodes.	Patients skin has cuts or abrasions in the area of electrode placement.	Discontinue use until skin is fully healed.
Skin becomes red by treatment.	Device has been used by the patient for too long periods.	Use device for a shorter treatment duration.
Treatment is uncomfortable, tingling sensations during use.	Intensity too high for required treatment.	Reduce intensity setting until treatment is comfortable. Note: Tingling sensations are normal and should be felt.

Programs Neurology

Progr.: 01		STIM + EMG alternating - Neurology				
		Phase 1				
Mode		STIM - EMG alternating				
Frequency Stim	Hz	45				
Pulse duration	μS	250				
Ramp up time	secs	1				
Ramp down time	secs	1				
Stim time	secs	4				
Rest time stim	secs	12				
Work time EMG	secs	8				
Rest time EMG	secs	12				
Number of cycles		15				
Progr.: 02		EMG - Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Mode		EMG				
Number of cycles		30				
Work time	secs	8				
Rest time	secs	12				
Progr.: 03		ETS - Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Mode		ETS				
Number of cycles		25				
Frequency work	Hz	45				
Pulse duration	μS	250				
Ramp up time	secs	1,5				
Ramp down time	secs	1,5				
Work time	secs	8				
Stim time	secs	4				
Rest time	secs	16				
Progr.: 04		Force 1 - Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	5		5		
Number of cycles			20			
Mode		W/R	EMG	Mod		
Frequency work	Hz	45		8-12		
Pulse duration	μS	250		300		
Modulation time	secs			10		
Ramp up time	secs	1				
Ramp down time	secs	1				
Work time	secs	6	6			
Rest time	secs	15	15			

Programs Neurology

Progr.: 05		Force 2 - Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	5		5		
Number of cycles			20			
Mode		W/R	ETS	Mod		
Frequency work	Hz	45	45	8-12		
Pulse duration	μS	250	250	250		
Modulation time	secs			10		
Ramp up time	secs	1	1			
Ramp down time	secs	1	1			
Work time	secs	8	8			
Stim time	secs		4			
Rest time	secs	16	16			
Progr.: 06		Force 3 - Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	5		5		5
Number of cycles			15		15	
Mode		W/R	EMG	W/R	EMG	Mod
Frequency work	Hz	45		45		8-12
Pulse duration	μS	250		250		250
Modulation time	secs					10
Ramp up time	secs	1		1		
Ramp down time	secs	1		1		
Work time	secs	6	6	6	6	
Rest time	secs	15	15	15	15	
Progr.: 07		Force 4 - Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	5			5	
Number of cycles			15	10		
Mode		W/R	EMG	ETS	Mod	
Frequency work	Hz	45		45	8-12	
Pulse duration	μS	250		250	250	
Modulation time	secs				10	
Ramp up time	secs	1		1		
Ramp down time	secs	1		1		
Work time	secs	6	6	8		
Stim time	secs			4		
Rest time	secs	15	15	16		

Programs Neurology

Progr.: 08		Relaxation - Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	5		5		
Number of cycles			20			
Mode		Mod	EMG	Mod		
Frequency work	Hz	8-12		8-12		
Pulse duration	μS	250		250		
Modulation time	secs	10		10		
Work time	secs		20			
Rest time	secs		10			
Progr.: 09		Rehabilitation 1 - Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	20				
Mode		W/R				
Frequency work	Hz	50				
Pulse duration	μS	300				
Ramp up time	secs	1				
Ramp down time	secs	1				
Work time	secs	10				
Rest time	secs	15				
Progr.: 10		Rehabilitation 2 - Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	5	15	10		
Mode		Cont	W/R	MOD		
Frequency work	Hz	5	60	8-12		
Pulse duration	μS	300	350	250		
Modulation time	secs			6		
Ramp up time	secs		2			
Ramp down time	secs		1,5			
Work time	secs		7			
Rest time	secs		14			

Programs Neurology

Progr.: 11	Strength Endurance + Explosive Force 1 - Neurology					
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	5	5			
Mode		W/R	W/R			
Frequency work	Hz	45	45			
Pulse duration	μS	250	250			
Ramp up time	secs	2	0,5			
Ramp down time	secs	1	0,5			
Work time	secs	5	2			
Rest time	secs	12	7			
Progr.: 12	Strength Endurance + Explosive Force 2 - Neurology					
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	5	5	5	5	
Mode		W/R	W/R	W/R	W/R	
Frequency work	Hz	45	45	45	45	
Pulse duration	μS	250	250	250	250	
Ramp up time	secs	2	0,5	2	0,5	
Ramp down time	secs	1	0,5	1	0,5	
Work time	secs	5	2	5	2	
Rest time	secs	12	7	12	7	
Progr.: 13	Partial peripheral paralysis - Neurology					
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	10				
Mode		W/R				
Frequency work	Hz	60				
Pulse duration	μS	450				
Ramp up time	secs	0,6				
Ramp down time	secs	1				
Work time	secs	4				
Rest time	secs	20				

Programs Neurology

Progr.: 14		Pain 1 - Urology + Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	30				
Mode		Mod				
Frequency work	Hz	70/110				
Pulse duration	μS	200/100				
Modulation time	secs	6				
Progr.: 15		Pain 2 - Urology + Neurology				
		Phase 1	Phase 2	Phase 3	Phase 4	Phase 5
Phase time	min	30				
Mode		Han				
100 Hz / 200 μS changing continuously with 2 Hz / 300 μS in 3:3-second rhythm						

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
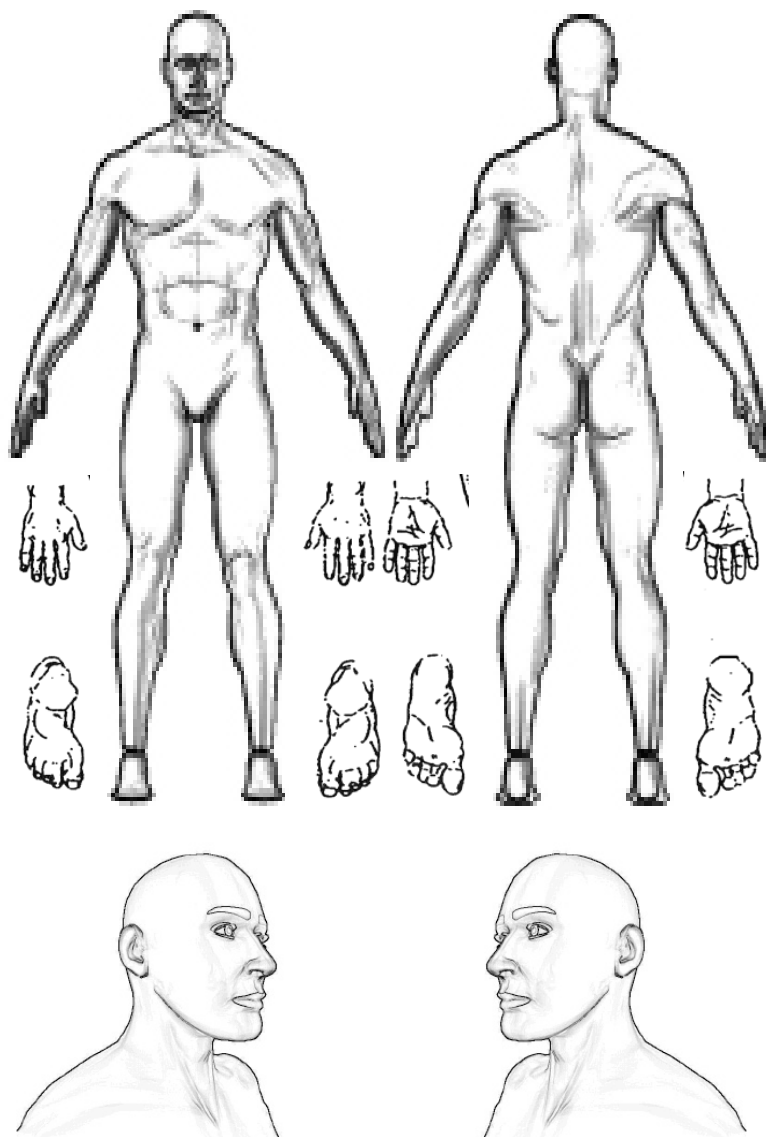
Please always use channel A first. If skin electrodes are used, both channels can be used.

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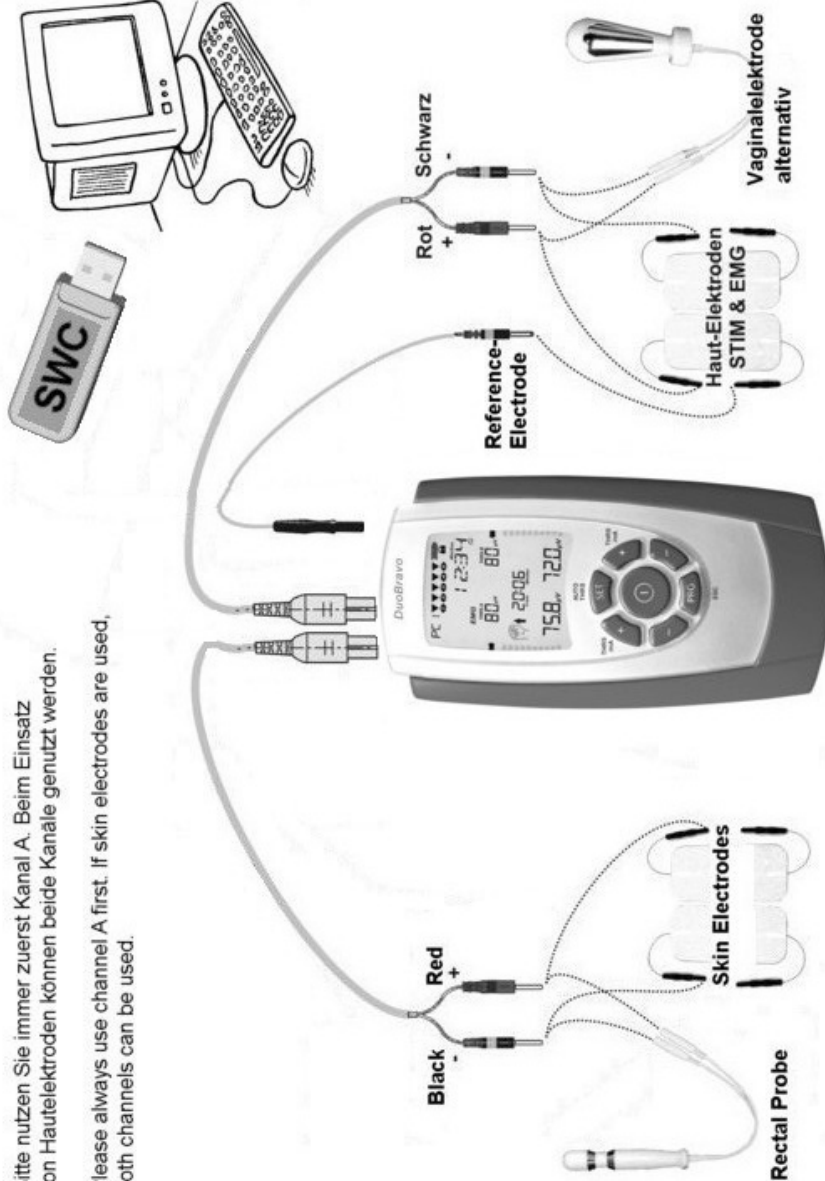In accordance with European regulations [1], used electrical and electronic equipment may not be discarded in household or unsorted municipal waste. Use separate collection facilities. Contact your local government or council for information regarding the collection systems available. If electrical appliances are disposed of in landfills or dumps, hazardous substances can leak into the ground water and get into the food chain damaging the health and well being of people and animals.

[1] European Community directive 2012/19/EU on waste electrical and electronic equipment.

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The batteries should be replaced when they are drained or stop functioning. The batteries supplied with the **DuoBravo** are recyclable. Remove the batteries from the unit and follow your local recycling guidelines.

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A list of literature regarding the topics of this manual can be found on our website „www.mtrplus.com“.

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Reserved for changes.

Printing or copying of this manual - also in an abridged version - is forbidden and can only be allowed with a written approval of “MTR+ Vertriebs GmbH”.

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Electromagnetic Compatibility

If the device is used according to the instructions given in this manual, it corresponds to the safety requirements for medical electrical systems IEC 60601-1. **DuoBravo** was checked and is in accordance with the limits set for medical products in the standard IEC 60601-1-2:2020. These limits serve to the safety of the user regarding electromagnetic disturbances. The **DuoBravo** is a Class B device and is suitable for home healthcare environments as well as typical clinical environments.

Guidance and manufacturer's declaration - electromagnetic emissions		
The DuoBravo unit is intended for use in the electromagnetic environment specified below. The customer or the user of the DuoBravo should ensure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The DuoBravo uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The DuoBravo is suitable for use in all establishments, including domestic and medical establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not Applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not Applicable	

Please adhere to these warnings when using the device:

Don't use **DuoBravo** near active HF surgical equipment or RF shielded rooms used for magnetic resonance imaging, where the intensity of EM disturbances is high.

Use of this equipment adjacent to, or placement on top of other equipment should be avoided in case output or EMG reading of the **DuoBravo** is effected. If this is necessary, Both the **DuoBravo** and the other equipment should be monitored to verify that they are both operating normally.

Use of accessories and lead wires other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used not closer than 4 meters to any part of the **DuoBravo** in EMG or ETS mode or 30 cm in Stimulation mode, including lead wires specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

Electromagnetic Compatibility

Guidance and manufacturer's declaration - electromagnetic Immunity			
The DuoBravo is intended for use in the electromagnetic environment specified below. The customer or the user of the DuoBravo should assure that it is used in such an environment			
Immunity Test	IEC 60601 Test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %
Electrical fast transient/burst IEC 61000-4-4	Power supply lines: ±2 kV input/output lines: ±1 kV	Not Applicable	DuoBravo is powered by batteries
Surge IEC 61000-4-5	line(s) to line(s): ±1 kV. line(s) to earth: ±2 kV. 100 kHz repetition frequency	Not Applicable	DuoBravo is powered by batteries
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° phases 0% U_T for 0.5 cycle And at single phase 0° 0% U_T for 1 cycle 70% U_T for 25/30 cycles 0% U_T for 250 cycles 0% U_T for 300 cycles	Not Applicable	DuoBravo is powered by batteries
Power frequency (50/60Hz) magnetic field IEC 61000-4-8	30 A/m 50 Hz / 60 Hz	30 A/m 50 Hz / 60 Hz	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Conducted RF IEC61000-4-6	150 KHz to 80 MHz: 3 Vrms 6 Vrms (in ISM and amateur radio bands) 80 % AM at 1 kHz	150 KHz to 80 MHz: 3 Vrms 6 Vrms (in ISM and amateur radio bands) 80 % AM at 1 kHz	Interference may occur in the vicinity of equipment marked with the following symbol: 
Radiated RF IEC61000-4-3	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz	
NOTE 1: U_T is the a.c. mains voltage prior to application of the test level. NOTE 2: At 80 MHz the higher frequency range applies. NOTE 3: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

Electromagnetic Compatibility


Guidance and manufacturer's declaration - electromagnetic Immunity

The **DuoBravo** is intended for use in the electromagnetic environment specified below. The customer or the user of the **DuoBravo** should assure that it is used in such an environment

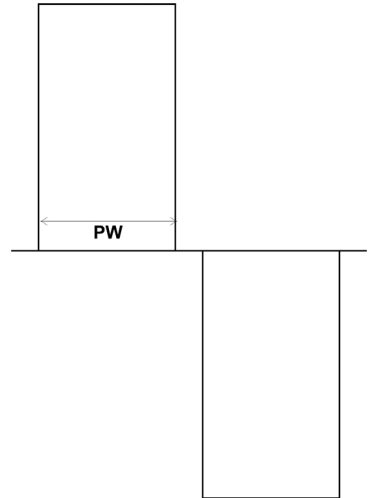	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Modulation (W)	Distance (m)	IMMUNITY TEST LEVEL (V/m)
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment)	385	380 – 390	TETRA 400	Pulse modulation 18 Hz	1,8	0,3	27
	450	380 – 390	GMRS 460, FRS 460	FM c) ± 5 kHz deviation 1 kHz sine	2	0,3	28
	710	704 – 787	LTE Band 13, 17	Pulse modulation 217 Hz	0,2	0,3	9
	745						
	780						
	810	800 – 960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation 18 Hz	2	0,3	28
	870						
	930						
	1720	1 700 – 1 990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation 217 Hz	2	0,3	28
	1845						
	1970						
	2450	2 400 – 2 570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0,3	28
	5240	5 100 – 5 800	WLAN 802.11 a/n	Pulse modulation 217 Hz	0,2	0,3	9
	5240						
	5785						

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